

# ScopeVAULT™


by ClinicalChoice

## Storage Cabinets



### ScopeVAULT™ ONE

- Colonoscopes
- Gastrosopes
- Duodenoscopes
- Bronchoscopes
- Enteroscopes

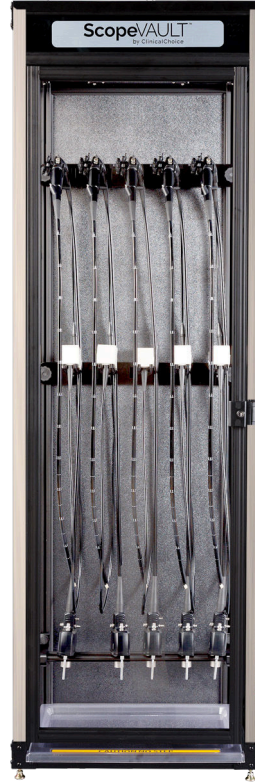

### ScopeVAULT™ PLUS

- ENT
- EUS (Ultrasound)
- EBUS (Ultrasound)
- TEE Probes
- Vaginal probes
- Rectal probes



Maximum Storage



Maximum Airflow



Small Footprint



Modular Design



Infection Control

## ScopeVAULT™ Cabinet Specifications

SIZE	WEIGHT	VOLTAGE
W 29", D 23", H 94"	250+/- pounds	120V 20A Circuit

ScopeVAULT™ storage cabinets and SafeCAP® endoscope transport and short-term storage system protect your valuable inventory and reduce cross-contamination.



Clinical Choice, LLC  
 204 Muirs Chapel Road  
 Suite 301Q  
 Greensboro, NC 27410

Phone: 877.572.9985  
 Fax: 877.572.9986  
[www.clinicalchoice.com](http://www.clinicalchoice.com)

## KEY FEATURES

- ✔ Universal design accommodates Olympus®, Pentax® and Fujinon® gastroscopes, colonoscopes, duodenoscopes, bronchoscopes, enteroscopes, EUS and EBUS (ultrasound) scopes. Specialty cabinets accommodate TEE, rectal and vaginal probes and other specialty endoscopes.
- ✔ 24/7 HEPA filtered air - 500 CFM airflow keeps your scopes clean and dry.
- ✔ No floor level vents to draw in dust and dirt.
- ✔ When the door is opened, bright LED lights illuminate your scope inventory to assist your staff.
- ✔ Top and bottom scope cradles secure your scopes as required by current regulations.
- ✔ Easy to clean drip tray can be removed with the door closed.
- ✔ Clear acrylic door enables you to visualize your entire scope inventory without opening the cabinet.
- ✔ Great use of space: small footprint, maximum storage.
- ✔ Modular design can be retrofitted to accommodate new and different scope inventory: e.g. ENT scopes and savory dilators.
- ✔ Earthquake brackets available.
- ✔ TJC inspection friendly. Meets or exceeds all current standards and guidelines including ANSI/AAMI ST91:2015, AORN, and SGNA.

## Order Information

PRODUCT ID	DESCRIPTION
SV-ONE	Standard GI cabinet upgradeable to channel drying.
SV-DRY	Standard cabinet with industry leading channel drying system installed.
SV-PLUS	Standard cabinet outfitted with purpose-built racks and cradles to accommodate specialty scopes (ENT, EUS, EBUS, TEE & PROBES). Can be outfitted with/upgraded to channel drying system.
SV-MOD	Our standard cabinet outfitted with custom racks and cradles to meet your exact storage needs. Can be outfitted with/upgraded to channel drying system.